**SPORT** **TEC**  
Physio & Fitness

==== [www.sport-tec.de](http://www.sport-tec.de) ====

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**Trimilin Trampolin Med Plus,  
ø 102 cm, bis 110 kg**

Art.Nr.: 22906

▶▶▶ [zum Produkt ...](#)

Trimilin Trampoline

▶▶▶ [zur Kategorie...](#)

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
Replace the springs as shown with a long screwdriver using the leverage effect

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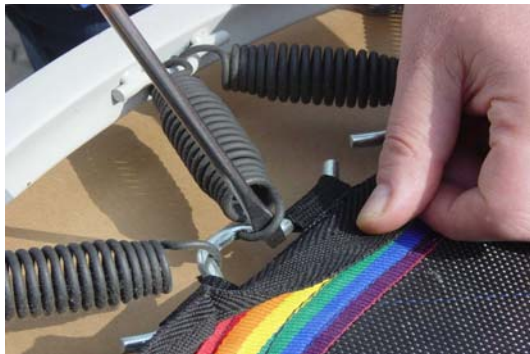
To dismount the springs use a screwdriver as shown. To mount the mat, place only 4 springs first, one each at 3, 6, 9 and 12 o'clock as shown. Then hook in the rest of the springs using the same technique with the screwdriver.

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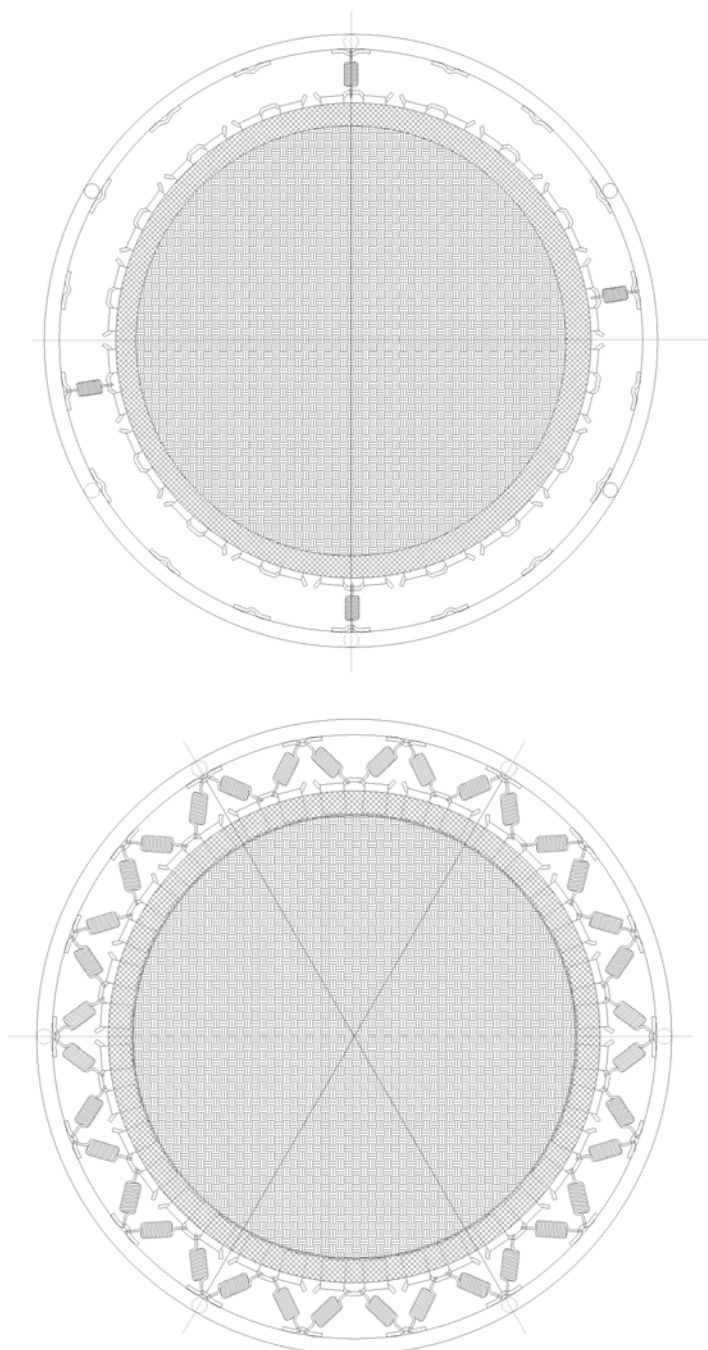## Mounting the frame cover of a Trimilin with foldable legs

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\*\*\*\*\*

Stand behind the trampoline which should be in a vertical position. The legs of the trampoline showing forward. Take the cover, the colour part outside, the backing showing to you. The legs remain folded.

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\*\*\*\*\*

Pull the rubber cable at the opening over the cap of the leg. Roll the trampoline sideways to the next leg and repeat pulling the rubber cable over the white cap. Do the same with all remaining legs. All of them stay folded

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\*\*\*\*\*

Carefully pull the cover up along the leg till the leg socket is fully visible. Then pull the cover over the frame to the other side.

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\*\*\*\*\*

Pull the cover along the other legs to the other side of the frame. On the other side of the frame the cover overlaps with 2-3 cm.

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\*\*\*\*\*

Carefully pull the cover completely over the frame of the trampoline. Repeat this with all remaining legs. Finally adjust the cover so that you have a regular circle and all the springs are covered.

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